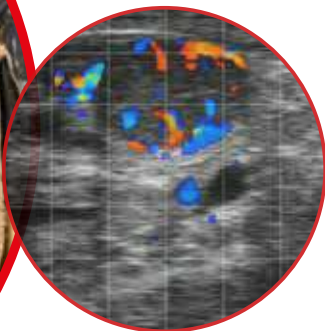




Easi-Scan:Go™

See in colour for better reproduction results.

The only wireless ultrasound scanner with colour flow doppler mapping software optimized for bovine reproduction.

CFM optimised for bovine reproduction

Colour flow doppler mapping (CFM) allows you to diagnose a negative PD from as early as 18 days post AI. Our CFM pre-set is designed with the bovine uterus in mind. Blood in the uterus does not circulate as fast as in other areas of the body and requires different imaging techniques, settings and scanner capabilities for accuracy. Easi-Scan:Go is the scanner that allows you to use CFM for bovine reproduction.

Improve KPI and output on the farm

CFM allows you to diagnose a negative PD from as early as 18 days post AI, therefore facilitating resynchronisation of the dairy cows up to 10 days earlier than before. This will save on average of £40 per cow and lead to greater milk yields. CFM also helps to determine the best recipient for embryo transfer. Getting the cow pregnant faster at the beginning of the breeding season will lead to greater calf weights in beef cows.

View the way you like it

Easi-Scan:Go is designed to make the scanning process easier and faster. You can view the image with the robust wireless colour BUGs, or via on the app on any mobile device. With the app you can also adjust CFM gain settings and view the image in both B and CFM mode simultaneously, making your scanning process more seamlessly productive.

facebook.com/IMVimaging
[twitter/IMVimaging](https://twitter.com/IMVimaging)
info@imv-imaging.com
 +44 (0) 1506 460023



**UP TO 3 SMART DEVICES
OR WIRELESS VGA BUGS**
can be connected to
view B mode, CFM or
CFM+B modes.

WATERPROOF
Sealed buttons with no
connectors

850 GRAMS

IMPROVED FLEXIBILITY
for better cable
management



USE CALIPERS
to measure the embryo
size and calculate foetal
age on the app.

RECORD VIDEOS
and save images
on your smart devices

**LOW COST LITHIUM-ION
BATTERY**
up to 5 h (B mode)
up to 4.5 h (CFM mode)

128 ELEMENTS
receiving on 8
channels, 26 frames
per seconds. Linear
and Curved probes
available.

Unique Technology of the Easi-Scan:Go.



At the forefront of bovine reproduction

We are making the latest technology accessible to all. Send your Easi-Scan:Go back to us for the CFM upgrade. Adopt the latest imaging technology, without learning how to use a new device!



Unrivalled Longevity

Tough farm conditions require reliable equipment, which is why we developed Shield-Tech, which makes the Easi-Scan:Go our most rugged scanner yet. Our fixed probe with strengthened cable has improved durability. The scanner can be fully immersed in a bucket of water improving biosecurity.



Superior image quality

Unique low power Smart-Beam Technology not only saves the battery life but also optimises the CFM image. The adapted low frame rate forms a superior image, which can be viewed in superb resolution on any wireless screen or colour viewing goggles.



Discover more at:

www.imv-imaging.com/easi-scan-go

www.sonorite.co.za



facebook.com/Sonorite



[twitter/IMVImaging](https://twitter.com/IMVImaging)



sonorite@mweb.co.za



+27 82 8880315